



## Gas Calibration Workflow Quick Guide

Vyntus™ ONE, Vyntus™ ONE-PFT, Vyntus™ ONE-DL, Vyntus™ BODY

1. Launch the **"Gas Calibration"** program - two ways to access the **"Gas Calibration"** program

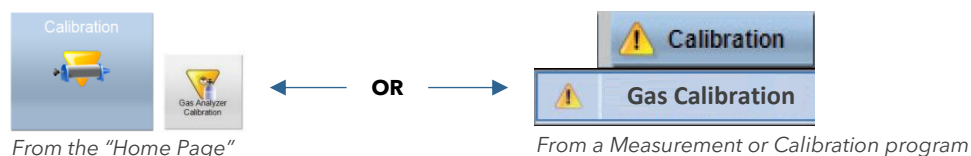

2. Connect the sample line to the Gas Calibration port on:

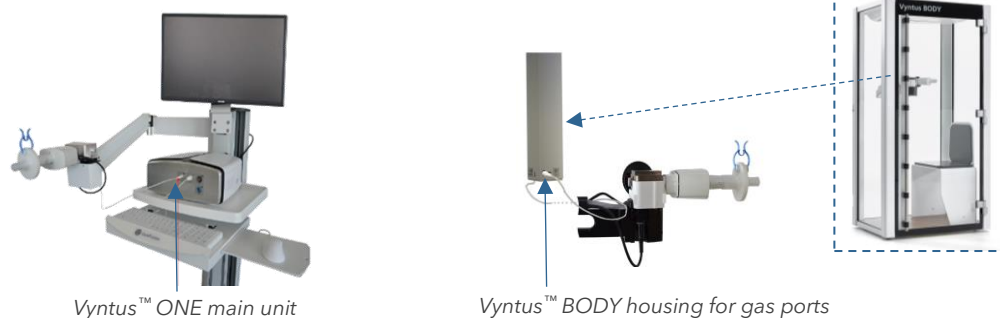

3. Click **<F1>** to **"Start"** the Gas Calibration



4. If the following window appears, select the option **"Calibration"**



5. After the Gas Calibration is completed, check if the results are within the limits = **"OK"**



6. Close the **"Gas Calibration"** program - two ways to close the **"Gas Calibration"** program

